

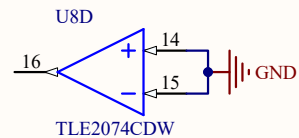
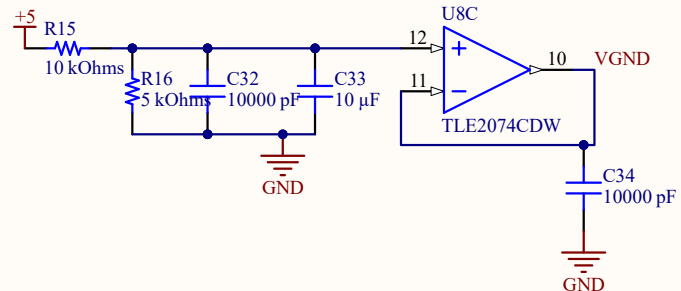
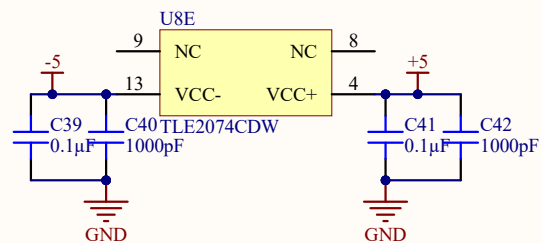


@VGND = 1.65 V,
so signal output is
DC biased to middle
of 3.3 V DR.

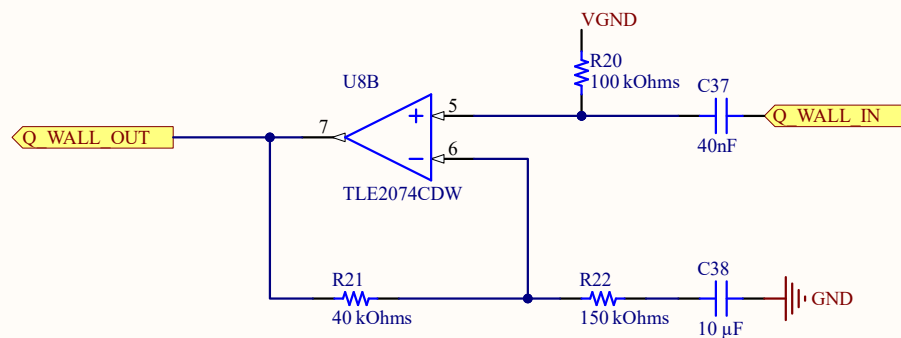

For ultrasound Doppler
measurements, the
arterial wall adds
low-frequency noise to
the spectral Doppler
waveform. The wall
filter shall remove these
low frequencies.

What cut-off frequencies
should I use?
22, 75, 128, 252 Hz
Component values give
 $f_c=27$ Hz

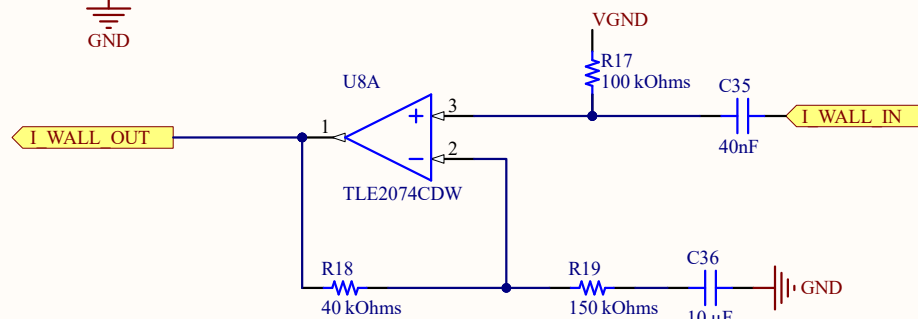

KAIST		(c) BMM Laboratory	
Project: ultrasoundsystem		Variant: [No Variations]	
Page Contents: 7_wallfilter.SchDoc		Checked by: -	
Size: A	DWG No: 0001		Revision: A
Date: 2023-05-03 JH		Sheet 6 of 6	